

A.M.I.® Soft Gastric Band System

Product Group Morbid Obesity | Issue 02/2010



> Setting Standards in Gastric Banding

A.M.I. Soft Gastric Band Premium

- Catheter is separate from safety closing mechanism
- Pre-curved shape reduces formation of folds
- Low pressure lowers risk of migration
- Three sizes accommodate full range of BMI's
- **NEW:** Ports with clip fixation

A.M.I.® Soft Gastric Band System

A.M.I. Soft Gastric Band Premium

The new generation of gastric bands is available in three different sizes, with a choice of four different port types!

Features:

- Gentle on enclosed tissue
- Low-pressure system
- Safety closing mechanism
- Unique catheter protection sleeve
- Catheter separate from closing mechanism
- 360° soft balloon

Advantages:


- Reduced risk of migration
- Wide adjustment range
- Easy closure
- No kinking, no leakage
- Less band slippage, less catheter damage
- Even distribution of pressure on stomach tissue

A.M.I. Soft Gastric Band System

Adjustable implant delivered with:

- Catheter
- Catheter protection sleeve
- A.M.I. Safety Port
- 2 x A.M.I. Port Needle
- Instructions for use
- Patient guide and patient card

1 implant, delivered sterile


	Standard Clip	Low Profile Clip	Standard	Low Profile
				
Height of port (without clip)	13 mm	9.8 mm	13.3 mm	9.8 mm
Short Band				
	AGB 372 CLIP	AGB 372_LP CLIP	AGB 372	AGB 372_LP
Standard Band				
	AGB 374 CLIP	AGB 374_LP CLIP	AGB 374	AGB 374_LP
Long Band				
	AGB 376 CLIP	AGB 376_LP CLIP	AGB 376	AGB 376_LP

A.M.I. Safe Puncture Ports



A.M.I. now offers the Soft Gastric Band Premium with a choice of four different ports. In addition to the conventional models in standard height or low-profile design, the bands are also available with revolutionary new clip-on ports in two different heights. Both these models feature a specially-designed base plate with an integrated clip, which can be slid under the fascia with ease for quick attachment of the port.

All the A.M.I. Safe Puncture Ports are made of bio- and MRI-compatible materials, and feature dual safety edges to prevent needle slippage. A special sleeve protects the catheter against puncturing and kinking, while the large silicone septum facilitates needle insertion.


NEW - Safe Puncture Port with Clip

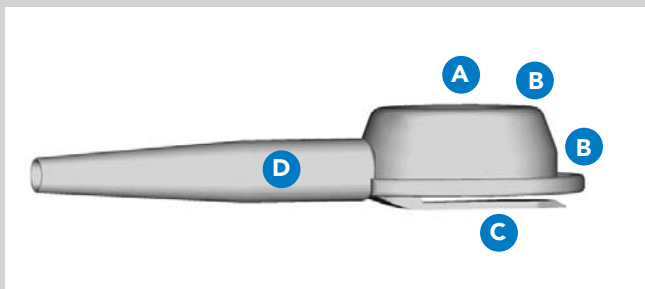

- A** Large silicone septum (12 mm)
- facilitates needle insertion
- B** Dual safety edges
- prevent needle slippage
- C** Clip
- saves sutures and operating time
- D** Catheter protection sleeve
- reduces kinking and leakage

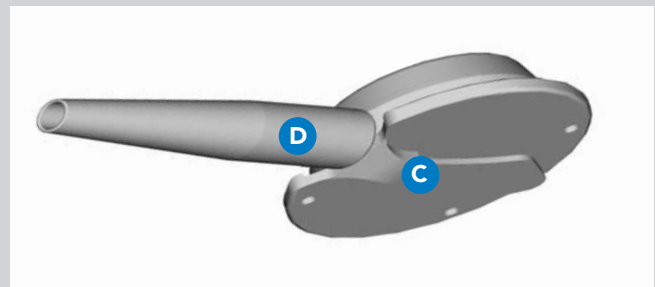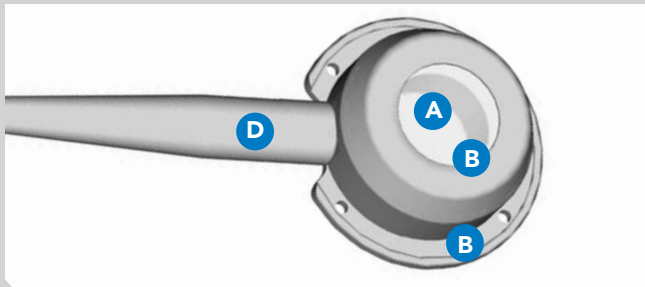

Years of Experience with A.M.I.'s Soft Gastric Band

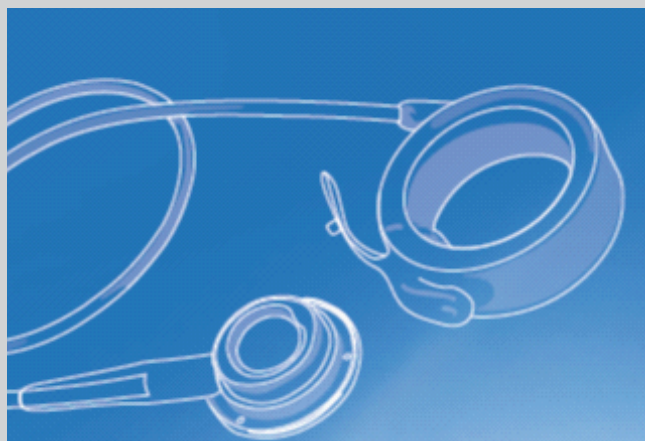








For more than 5 years we have been producing the A.M.I. range of "premium" gastric bands.

As with many other products, we set a trend with these bands which proved to be years ahead of their time. In particular their softness set our bands apart from those of our competitors; the same softness which is still being gentle on patients today.

And that softness is also the reason you and your patients can look forward to:

- extremely low erosion rates
- extremely low migration rates, and
- an extremely low incidence of leakage!

A.M.I.® Morbid Obesity Accessories

Order Code	Morbid Obesity Accessories	Technical Details
AGB 355 	A.M.I. Gastric Tube Silicone gastric tube with eccentric balloon to measure the proximal gastric pouch	1 tube, delivered sterile 39 French Length: 735 mm Inner diameter: 8 mm Outer diameter: 13 mm
AGB 379 	A.M.I. Gastric Band Tubing Connector Titanium connector piece for repositioning of port or catheter Delivery includes spare piece of catheter	Delivered non-sterile, steam autoclavable Length: 25 mm
AGB 382 	A.M.I. Port Needle Special needle with lateral hole. For safe, non-cutting punctures of the port septum	25 needles / box, delivered sterile Size: 20 G Length: 50 mm
ASP 310 	A.M.I. Safe Puncture Port Standard-height port of titanium, large base plate, safety edges and catheter protection sleeve	1 port, delivered sterile
ASP 300 	A.M.I. Safe Puncture Port Low Profile Low-profile port of titanium, large base plate, safety edges and catheter protection sleeve	1 port, delivered sterile
ASP3021 	A.M.I. Safe Puncture Port Clip Standard-height port of titanium, large base plate, safety edges, catheter protection sleeve and integrated clip	1 port, delivered sterile
ASP3011 	A.M.I. Safe Puncture Port LP Clip Low-profile port of titanium, large base plate, safety edges, catheter protection sleeve and integrated clip	1 port, delivered sterile

A.M.I. Headquarters:

A.M.I. GmbH
 Im Letten 1
 6800 Feldkirch . Austria
 t +43 5522 90505-0
 f +43 5522 90505-4006
 e info@ami.at
www.ami.at